

Confidence in your water quality analysis. Anytime. Anywhere.

Hach®'s Mobile Sensor Management provides clarity through easy access to crucial information, helping you manage your process and be proactive in your maintenance.



Guidance through actionable sensor information




Focus on your highest priorities, avoid unexpected downtime and reduce maintenance costs with access to actionable information and proactive notifications about maintenance needs.

Alignment of process and laboratory measurements

Verify and adjust process instrument measurements using laboratory reference values with the tap of a finger.

Simple maintenance instructions in the palm of your hand

Perform instrument maintenance quickly and accurately by following step-by-step instructions on your mobile device.

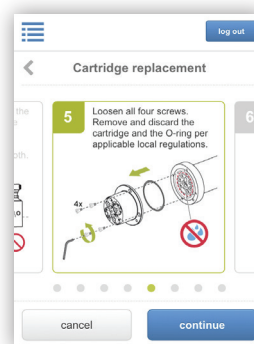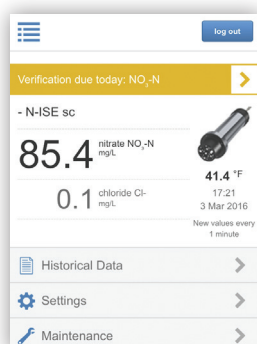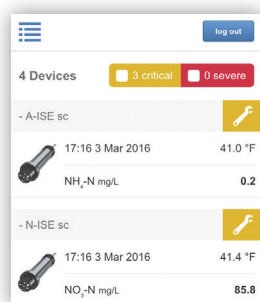


Principle of Operation

Mobile Sensor Management web-based software provides access to the enrolled Hach Analytical Instrumentation (1) through any Internet enabled device ("Peripheral Device") with a standard Internet Browser (2).

All access to the Mobile Sensor Management Web Server (3) is controlled by secure authentication. The status information provided gives an overview of the performance of each individual analytical instrument. It also provides specific information about upcoming maintenance activities for each instrument.

Maintenance procedures can be started easily from the dashboard of the Mobile Sensor Management software. Interactive and step-by-step guidance through the various maintenance steps allows for quick and accurate execution. Integrating Hach's DR3900/DR6000 Spectrophotometer (4) enables simple and convenient verification and calibration of the Analytical Instrumentation (1).

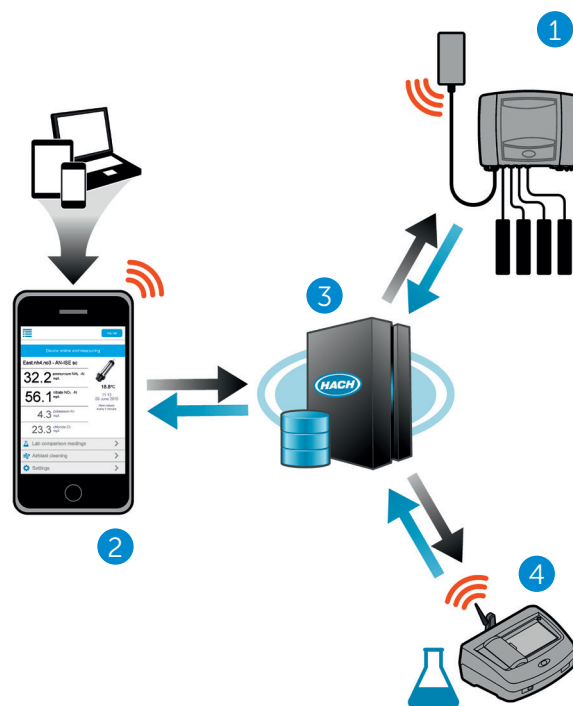

The communication between the Hach instruments (1) and (4) is established through Wi-Fi or a 3G/4G modem.

System Components

Compatible instruments:

- | | | | | | |
|-------------|------------|-------------|---------------|------------------|-------------------|
| • 1200-S sc | • 9184sc | • AN-ISE sc | • LDO 2 sc | • Phosphax sc LR | • SS7 sc |
| • 3400sc | • 9185sc | • CL10sc | • N-ISE sc | • SC1500 | • TSS sc |
| • 3700sc | • 9187sc | • CL17sc | • Nitratax sc | • SC4200c | • TU5300/TU5400sc |
| • 3798-S sc | • A-ISE sc | • DR3900 | • pHd sc | • Solitax sc | • Ultraturb sc |
| • 8362sc | • Amtax sc | • DR6000 | • Phosphax sc | • Sonatax sc | • Uvas sc |

Be confident with Hach Service

With Hach Service, you have a global partner who understands your needs and cares about delivering timely, high-quality service you can trust. Our Service Team brings unique expertise to help you maximise instrument uptime, ensure data integrity, maintain operational stability, and reduce compliance risk.